

Syringe Pump Communicator

User's Manual

Installation

Insert CD ROM, a file containing the CD's contents will open automatically. Right click on "Syringe Pump Communicator". Select "Send To" then "Desktop" to save the application on your desktop.


*Compatible with Windows 2000, XP, Vista, & Windows 7

1. Connect the pump to the PC using either an RS-232 or USB cable
2. Turn pump on
3. Start Communicator Program
4. Select the proper COM port.
5. Select the Baud Rate that matches the pump's Baud Rate (RS-232 only)

Overview

Pump Communicator allows the user to:

- Upload programs from the pump to a PC
- Download programs from PC to Pump
- Send commands to the pump using a PC

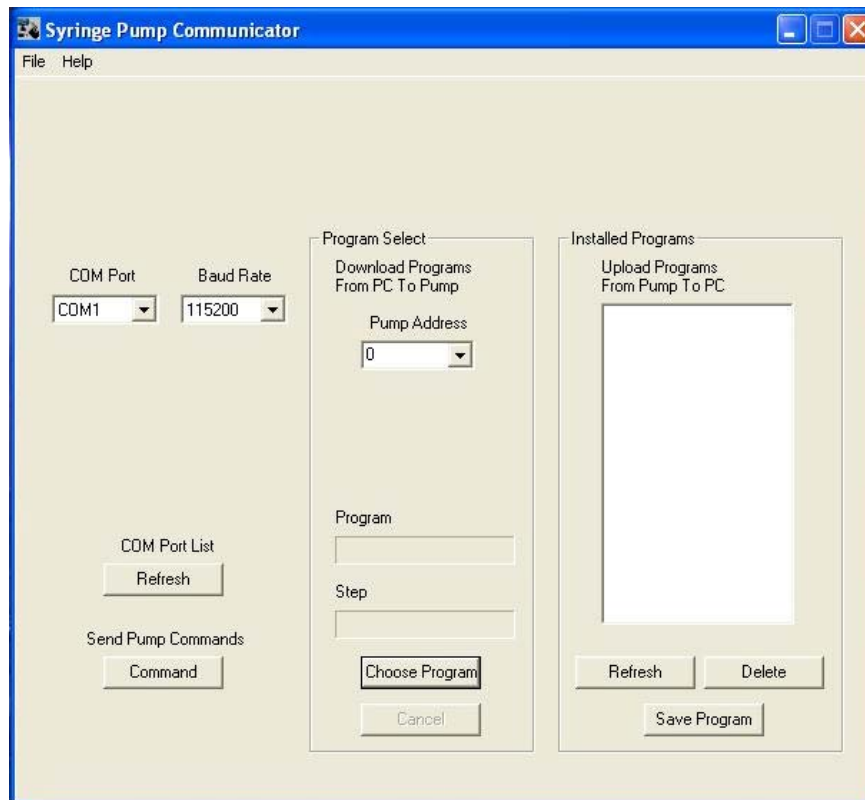

Pump Command *
Line

PC → Pump
Program Control
(programmable
pumps only)


Pump → PC
Program Control
(programmable
pumps only)

*No Need for HyperTerminal

Upload – Pump to PC

1. Select Program(s)


2. Click 'Save Program'

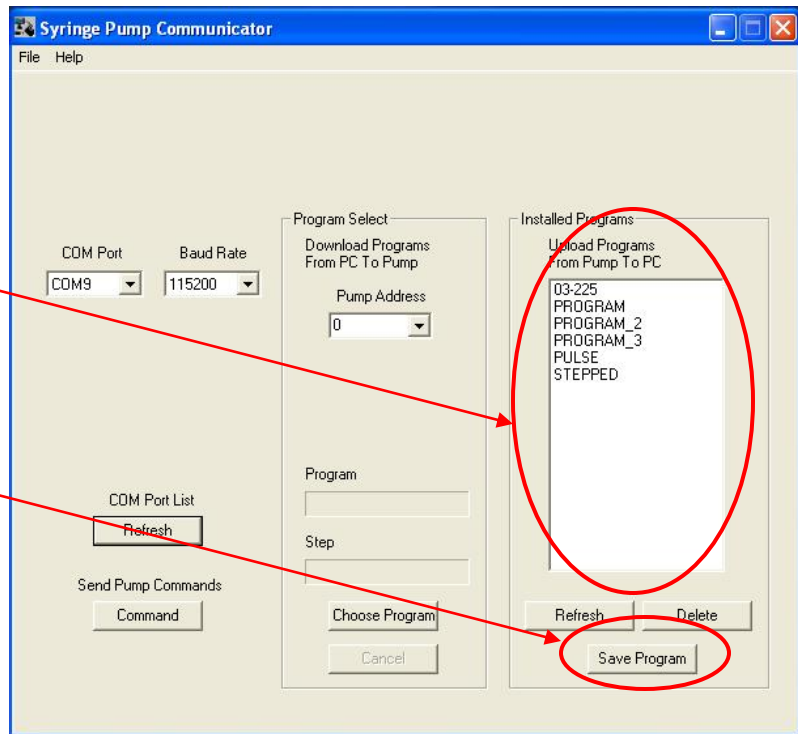

NOTE: Hit 'Refresh' to update the list after creating or editing programs on the pump

3. Select the drive and folder location

4. Click 'Save'

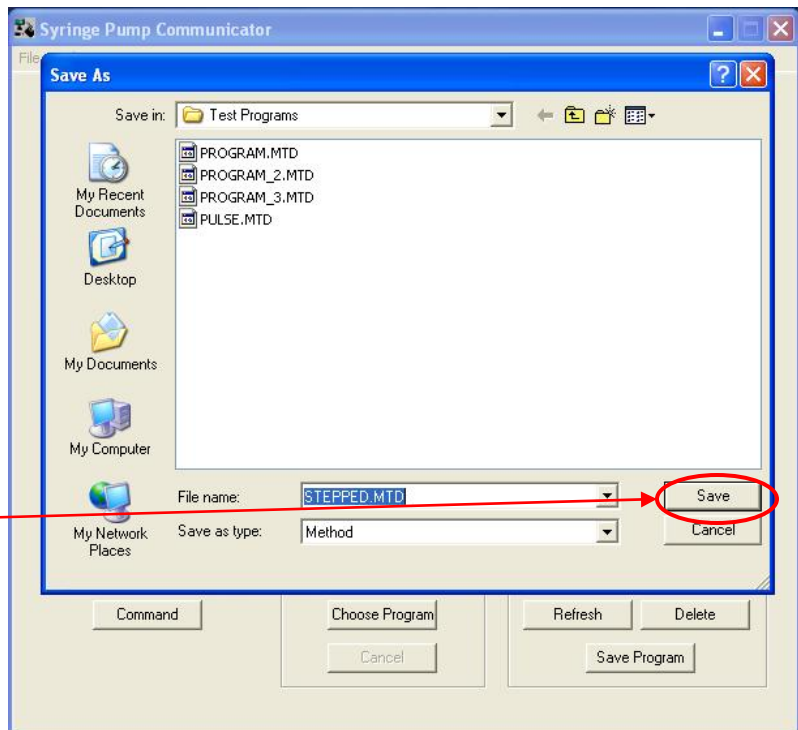


*Do not upload/download programs while pump is running
Communicator Program works best while pump is idle on Configuration screen


Download – PC to Pump

1. Click 'Choose Program'

2. Select Program(s) to download

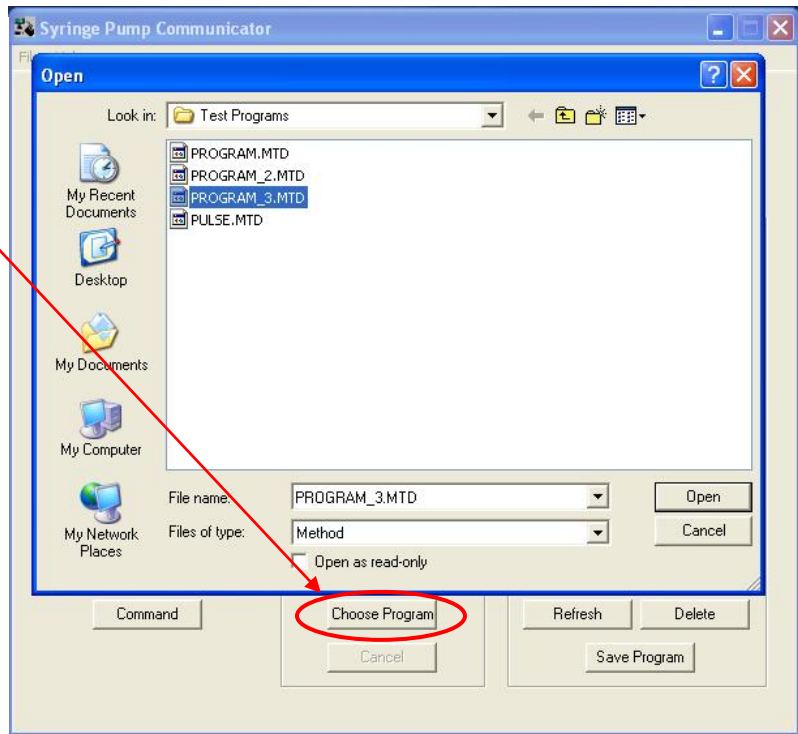


3. Program name and step number being downloaded is shown

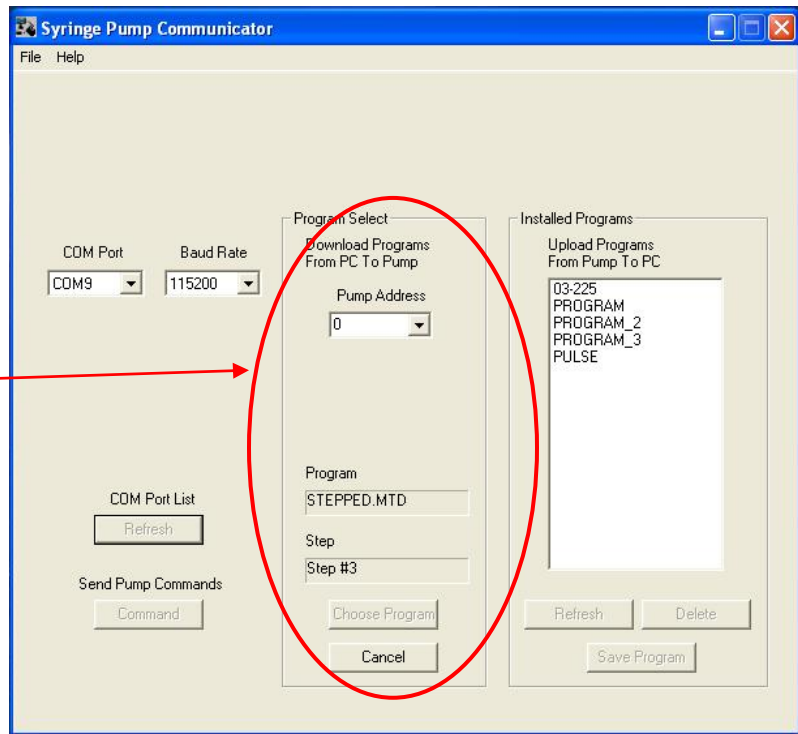

*Do not upload/download programs while pump is running
Communicator Program works best while pump is idle on the Configuration screen

If a program already exists on the pump, Syringe Pump Communicator will not overwrite the file if you attempt to re-download it

NOTE: Changing a file name, on the PC, will not change the program name within the file

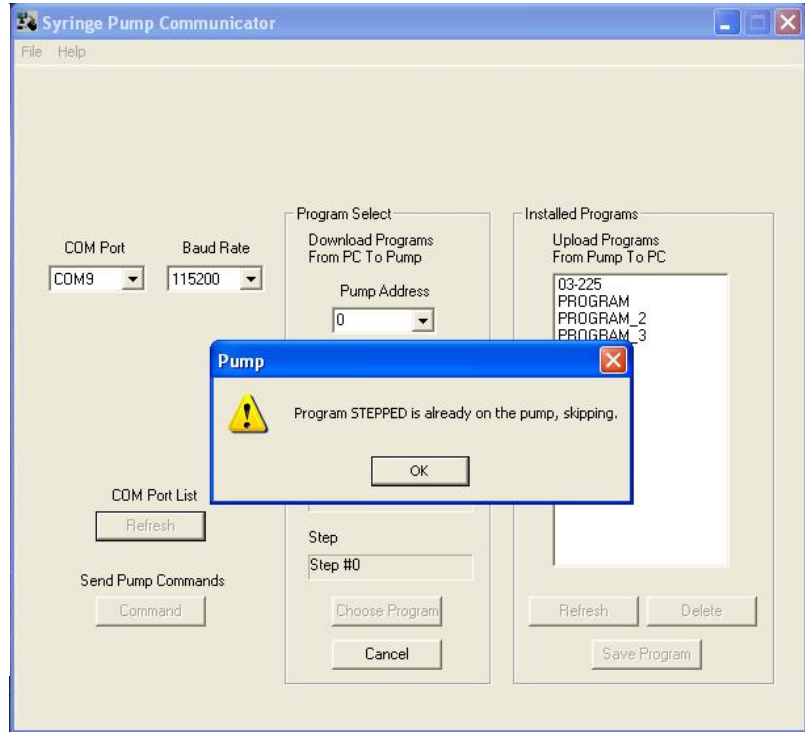


Sending Commands to the Pump

1. Click 'Command'
2. Type commands
3. Click 'Send'

Note 1: See pump's user manual for complete list of commands

Note 2: If unexpected responses are seen type 'echo off'

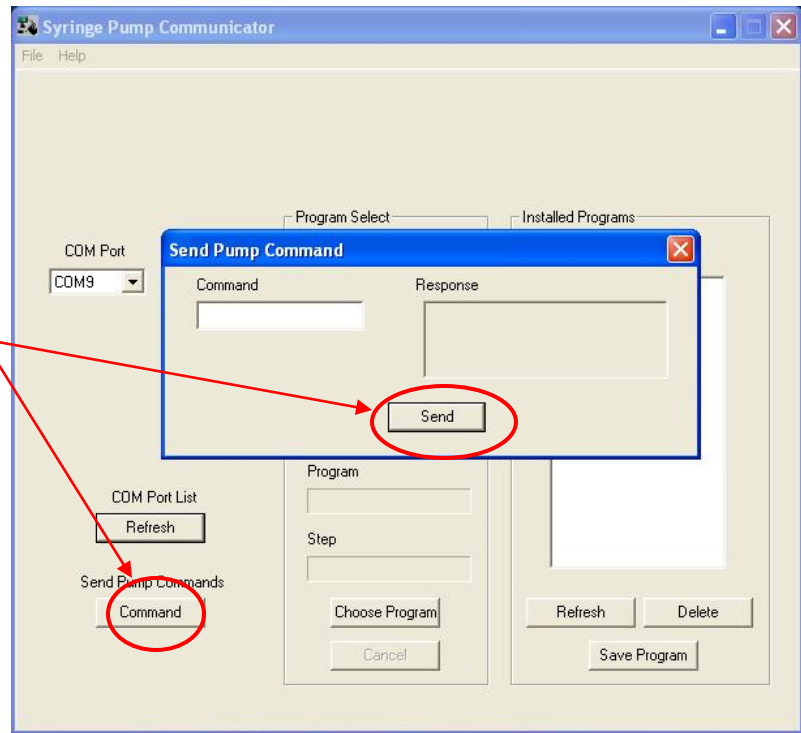

Troubleshooting

1. Do not upload/download programs while pump is running
 - Communicator program works best while pump is idle on the Configuration screen
2. Order of Operations:
 - Connect USB or RS-232 cable
 - Turn pump on
 - Start Syringe Pump Communicator program
3. If using RS-232, check Baud rate – pump vs. program
4. Message 'Pump not Programmable'
 - If not a programmable pump – contact pump manufacturer for an upgrade
 - If pump is programmable – go to command line and type 'echo', if echo is on type 'echo off'